

BEFORE & AFTER



BEFORE



AFTER




BEFORE




AFTER




BEFORE

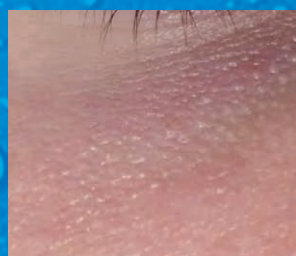

AFTER

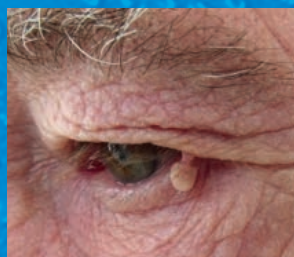

BEFORE

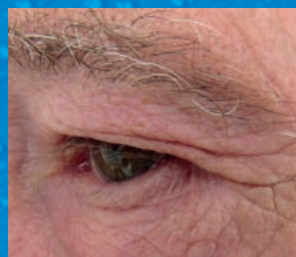

AFTER

CryoPen®

JOIN THE COLD
REVOLUTION

BEAUTY BY THE COLD™



MD *Aesthetics*
MEDICAL

CRYO  STHETICS

www.cryosthetics.co.uk

CRYO  STHETICS

WHAT IS CRYOTHERAPY AND WHAT DOES IT TREAT?

CryoPen™ is an advanced cryotherapy innovation that is a fast, effective, safe and new solution for removal of skin imperfections. Procedures typically last from 5 to 30 seconds, so you can be in and out of the salon quickly.

The CryoPen™ emits a fine jet of nitrous oxide under high pressure, which allows the clinician to work with millimetre precision. This destroys the tissue by freezing the inter-cellular fluid, forming ice shards and crystals which rupture the membrane, thereby destroying the cell. That means there will be no collateral damage to healthy tissue. It's so incredibly accurate and practical.

Cosmetic procedures such as skin tags, milia and cherry angiomas are no longer being treated by many doctors on the NHS. With the CryoPen™ you can now have these treatments in the beauty salon.

CryoPen™ offers the ideal treatment for:

- Skin tags
- Solar lentigo
- Age spots
- Milia
- Warts
- Cherry angioma
- Viral verrucae



HOW IT WORKS

The applicator is held as close as possible to the skin imperfection and moved quickly towards and away from it. This could be from 1-30 seconds, depending on the size and depth.

- Removes unsightly skin tags, warts and age spots for that perfect finish
- Minimally invasive



FREQUENTLY ASKED QUESTIONS

Are there any side effects?

The area treated is frozen by cryotherapy. This kills the localised cells in the area which cause them to shrivel and die. The treated area may go darker over a few weeks, but will return to normal pigmented skin in about 6 weeks.

How long are CryoPen™ treatments?

The duration of the treatment time will depend on the surface area and the thickness of the tissue being treated. Upon the clinical evaluation of the depth of the area, the clinician will choose an application time from 1 to 30 seconds.

How much will the treatment cost?

The cost of treatment will depend on the clinic and the size of the area. Make an appointment with your clinician to discuss costs and to assess your area.

How safe is cryotherapy?

Cryotherapy is a relatively low-risk procedure. Cryopen™ delivers N2O directly to the area, and not the healthy surrounding tissue. CryoPen™ is both CE and FDA registered.

Is the treatment painful?

There will be a pain sensation similar to a stinging nettle on the skin when the nitrous oxide reaches the bottom of the area. There might be a little residual stinging for a few minutes after treatment.